

Enhancing Diagnostic Confidence: Unleashing the Power of Iodine Mapping for CTPA Imaging

"We use iodine mapping for CTPA imaging routinely for enhanced detection of pulmonary embolism. It is particularly valuable in identifying sub-segmental emboli as well as chronic thromboembolic disease. The additional information provided by iodine maps improves diagnostic accuracy and aids in making crucial treatment decisions."

Dr. Russell Bull


Royal Bournemouth Hospital, United Kingdom



Patient History

Patient arrives to the emergency department with shortness of breath. A CT pulmonary angiogram (CTPA) is performed. The iodine maps are automatically generated and sent to PACS for immediate review.


Result

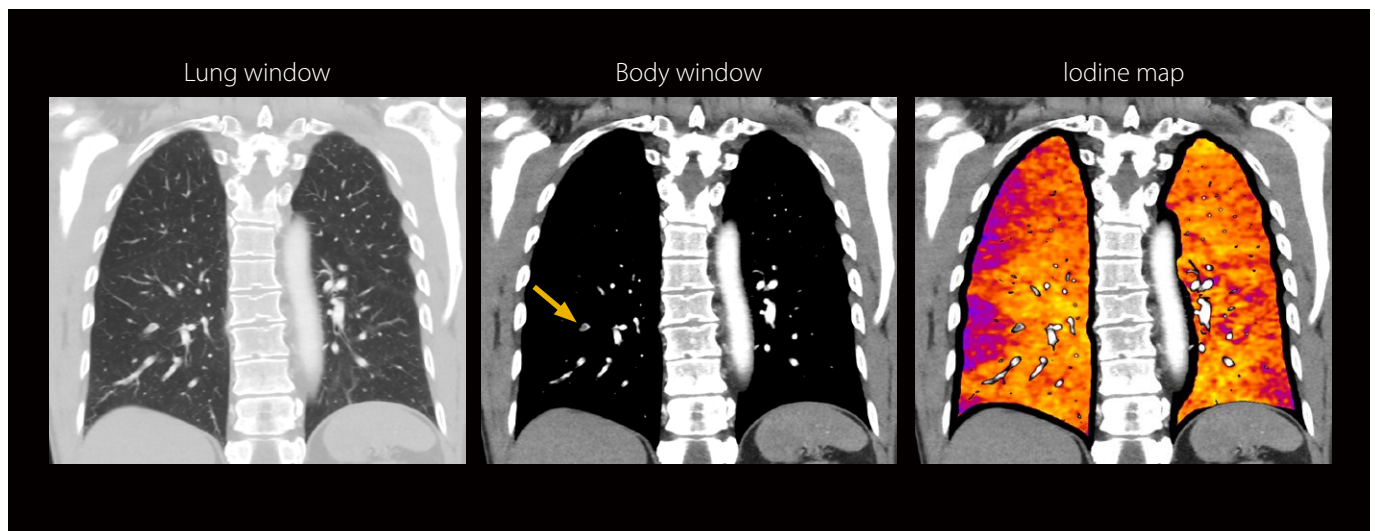

The iodine maps demonstrate an area of hypo-perfusion indicative of a pulmonary embolus. This allowed the occlusive embolus to be rapidly identified in the right pulmonary circulation on the conventional gray scale anatomical images. According to Dr. Bull, the added information from the iodine maps facilitates the expedited and accurate diagnosis of pulmonary embolism, potentially saving valuable time in patient care.


Technology

Iodine mapping SCT utilizes a 3D deformable registration algorithm for precise iodine signal extraction, compensating for patient motion. The extracted iodine is displayed in color, fused with the contrast-enhanced image, and automatically sent to PACS.

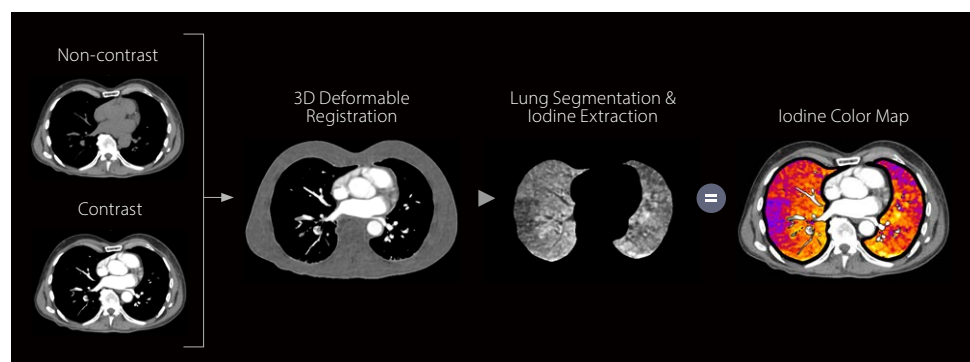

Subtraction Viewer

When using the SCT Viewer* it is also possible to interactively scroll the color iodine maps synchronized with the grayscale images for quick visualization and detection of hypoperfused areas.

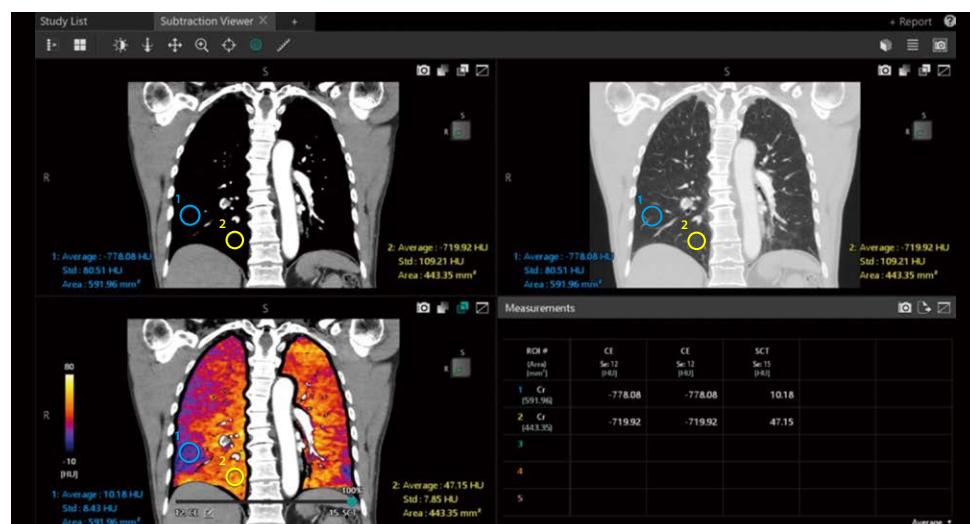

Conclusion

In this patient the iodine maps demonstrated an area of hypo-perfusion allowing the rapid identification of a pulmonary embolus. Iodine mapping boosts diagnostic confidence, guiding informed decisions for pulmonary vascular disease management.

All of Canon Medical's CT scanners can provide this unique package, enabling color blood flow maps for any single source, single energy, multiphase protocol with a pre-contrast scan.

*Subtraction is an option on the CT console and the Viewer is an option on the Vitrea Workstation

CANON MEDICAL SYSTEMS CORPORATION

<https://global.medical.canon>

©Canon Medical Systems Corporation 2024. All rights reserved.
Design and specifications are subject to change without notice.

Model number: TSX-307A MCSCT0011EAB 2024-08 CMSC/Produced in Japan

Canon Medical Systems Corporation meets internationally recognized standards for Quality Management System ISO 9001, ISO 13485. Canon Medical Systems Corporation meets the Environmental Management System standard ISO 14001.

Clinical results may vary due to clinical settings, patient preparation, and other factors. The views and opinions expressed in this brochure are those of the clinicians and do not necessarily reflect the views of Canon Medical Systems.

Acquisition

Scanner Model: Aquilion Serve
Scan Mode: Ultra Helical Subtraction
Collimation: 0.5 mm x 80
Exposure: 120 kV, ^{SURE}Exposure
Rotation Time: 0.35 s
CTDIvol: 1.6, 2.0 mGy
Total DLP: 126.9 mGy-cm
Effective Dose: 1.78 mSv
k-factor: 0.014*

Key Benefits

- Available with any single source, single energy CT
- Automatic image processing and transfer
- Display local iodine distribution
- Detect/rule out contrast enhancement
- Detect local perfusion differences

Key Benefits

- Automatic hanging layouts
- Synchronized review of all displayed images
- Geometric and HU measurements

* American Association of Physicists in Medicine (AAPM) Report 96, 2008.

Made For life